

HPLC Data Visualizer

v1.3.0 Quick Guide

↓ Download SVG figure



Windows portable edition

1. Start and import

Windows portable edition

- Fully extract the ZIP, then double-click HPLC Data Visualizer.exe. Python is not required.
- If SmartScreen appears on first launch, choose More info, verify the file source, and select Run anyway.
- Closing the final browser page exits the app after about 10 seconds. You can also use Close local app in the sidebar.

Supported data

- CSV: time in the first column and signal in the second; comma, semicolon, and tab delimiters; UTF-8, UTF-16, and GB18030.
- A CSV may contain several sequential two-column signals. Intensity is shown as likely UV/DAD and Counts as likely MS TIC/EIC.
- CTX: time;signal; numeric data from the first line, or data inside a [Chromatogram Data] section.

Multi-signal CSV

The three preset buttons only change the checkboxes. Review the selection and click Confirm and plot to import. Detector labels are inferred from column names and are provided as hints.

Choose detector signals to import

Detected 1 files containing multiple signals. Select the signals to plot first; detector types are inferred from column names only.

Import all

UV/DAD only

MS only

1.csv

☒ Intensity — Likely UV/DAD signal

☒ Y(Counts) — Likely MS TIC/EIC signal

Confirm and plot

2. Organize curves and set figure width

Order, names, and diagnostics

- Drag sample cards up or down to reorder curves.
- Expand Rename or remove curves to edit display names used by legends, curve labels, and SVG exports.
- Expand Import diagnostics to review columns, signal type, point count, time range, skipped rows, duplicate times, encoding, and delimiter.

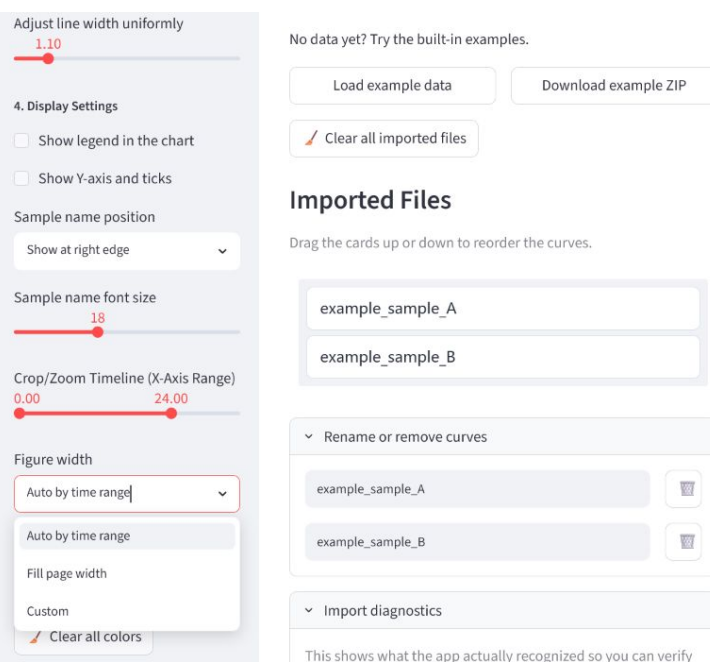


Figure width

- Auto by time range sizes the output from the current X-axis span and keeps short ranges compact.
- Fill page width follows the available browser width and exports at the current rendered size.
- Custom uses an explicit pixel width. A minimum width is retained for the chart and sample labels.

3. Chart interaction and export

Start with direct mouse controls

Use the mouse wheel to zoom, drag the chart to pan, click a curve to add or remove a retention-time marker, and double-click to restore the initial zoom. These actions do not require a toolbar button first.

Gray toolbar buttons

	Zoom Drag a rectangle for a precise zoom. Mouse-wheel zoom works without selecting this button.
	Pan Return to normal click-and-drag panning after using a selection tool.
	Box Select Drag across a peak region to add peak-area highlighting.
	Reset axes Restore the range set in the sidebar. Double-clicking the chart also restores the zoom.
	Download plot Download the current chart as SVG. The large button above the chart also exports SVG.

Peak highlighting

- Select Box Select, then drag across the peak region to highlight it.
- Select Pan when finished to return to normal dragging. Highlight targets can be changed or cleared from the sidebar.

Export

- The large Download SVG figure button and the toolbar download button both export SVG.
- Before exporting, check the time range, figure width, sample-name position, and font size.